



CATHETERISATION TRAINERS

60150

Advanced, Complete Set

60151

Advanced Male

60155

Advanced Female

60166

Basic, Value Set



Learning has never been more lifelike


limbsandthings.com

ADVANCED CATHETERISATION TRAINER - COMPLETE SET

60150 Complete Set

Simulated Patient

Our Advanced Catheterisation Trainer facilitates the teaching of urethral and supra-pubic catheterisation, as well as how to demonstrate self-catheterisation to patients (when using the optional Self Catheterisation Stand). It has interchangeable accurate male and female anatomy, and both sets of genitalia are realistically soft and pliable. The feel of the catheter passing along the urethra into the bladder corresponds closely to real life.

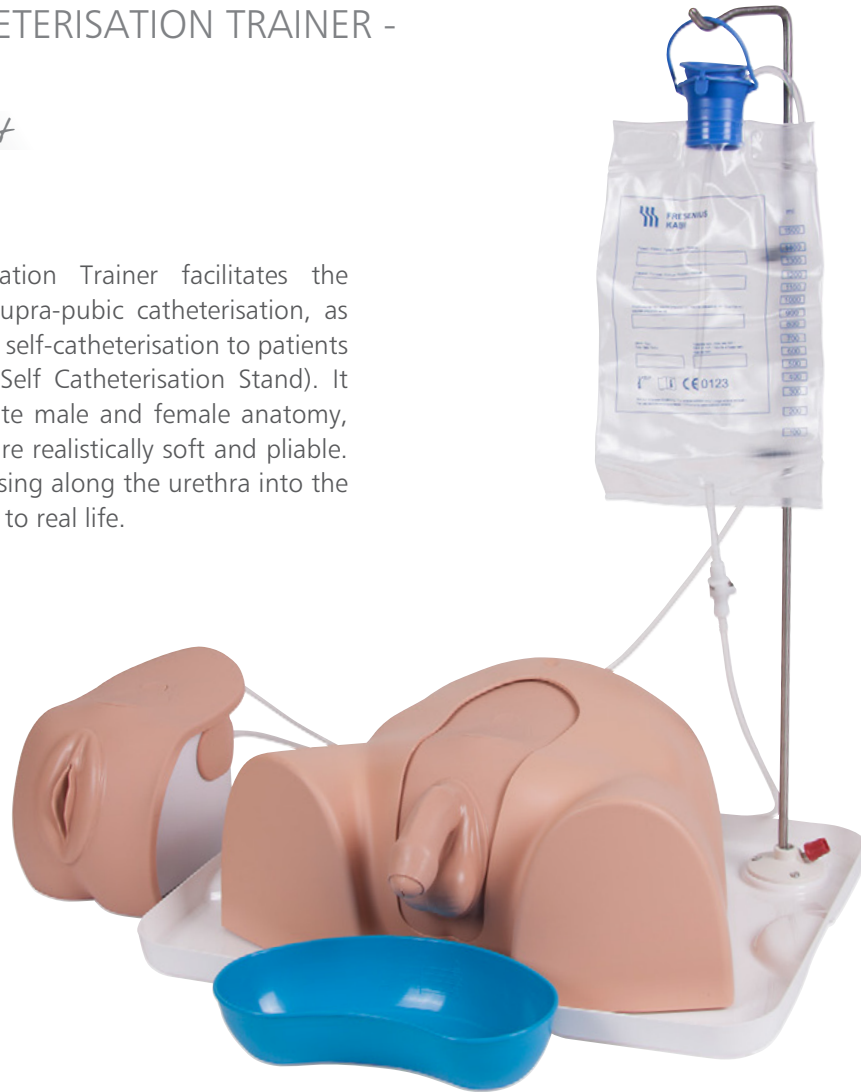

CATHETERISATION TRAINERS

Skills Gained

- Correct handling of male and female anatomy
- Aseptic catheterisation technique
- Catheter placement: accepts catheters from 14-16 French (male) and 12-16 French (female)
- Fluid management
- Withdrawal of catheter
- Intermittent Self Catheterisation (ISC) when used with the optional Self Catheterisation Stand (**60167**)
- Suprapubic catheter insertion and catheter management using Suprapubic Bung*

* Advanced Trainers only

Product Qualities

Anatomy

- Male anatomy: realistic meatus, flaccid penis and replaceable foreskin
- Female anatomy: partable labia with realistic urethral meatus and vaginal opening

Realism

- Supple urethra and resistant sphincter providing realistic response

Key Features

- 1 litre fluid bag on stand provides pressurised fluid flow
- Syringe supplied with water-based lubricant to simulate proprietary local anaesthetic gel
- Non-drip valve
- Reusable double-sleeve catheter packaging is supplied for teaching aseptic catheter handling

Versatility

- Easily interchangeable male and female catheterisation units

Safety

- Product is latex free (catheters provided contain latex)

Cleaning & Care


- Skin surface is washable using soap and water

ADVANCED CATHETERISATION TRAINERS - MALE & FEMALE



The Advanced Male and Female Catheterisation Trainers are available separately. Both are supplied with a suprapubic site and bung.

BASIC CATHETERISATION TRAINERS - VALUE SET

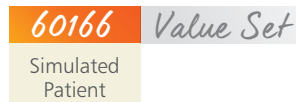

A cost effective Value Set containing basic versions of each of the Male and Female Catheterisation Trainers. These do not contain the suprapubic bungs.



Contents

	60150	60151	60155	60166
• 60152 Male Catheterisation Unit	✓	✓	✗	✓
• 60154 Catheterisation Foreskin (x3)	✓	✓	✗	✓
• 60156 Female Catheterisation Unit	✓	✗	✓	✓
• 60158 Suprapubic bung (x2)	✓	✓	✓	✗
• 60161 Standard bung	✓ (x2)	✓	✓	✓ (x2)
• 60159 Male Catheterisation Sphincter Kit	✓	✓	✗	✓
• 60162 Female Catheterisation Sphincter Kit	✓	✗	✓	✓
• 60160 Aseptic Catheterisation Sleeves (x2)	✓	✓	✓	✓ (x2)
• 50183 Mock LA Gel with Syringe	✓	✓	✓	✓ (x2)
• 50185 Antifungal Agent (Milton Sterilising Fluid)	✓	✓	✓	✓ (x2)
• Drip tray with bag stand	✓	✓	✓	✓ (x2)
• 1 Litre fluid bag & plastic kidney dish	✓	✓	✓	✓ (x2)
• 14 French foley catheter (male)	✓	✓	✗	✓
• 12 French foley catheter (female)	✓	✗	✓	✓
• Pelvic shell	✓	✓	✓	✓ (x2)
• 2 Litre jug	✓	✓	✓	✓ (x2)
• Carry case	✓	✓	✓	✓ (x2)



Accurate anatomical modelling



Aseptic catheterisation technique



Easily interchangeable units

Optional Extras & Consumables

SUPRAPUBIC CATHETERISATION MANAGEMENT BUNG

60168

This Pre-Ported Suprapubic Bung allows for the management of catheters after surgical procedure has taken place.



SELF CATHETERISATION STAND

60167

Use with Catheterisation Trainers: 60150, 60151, 60155, 60166 for teaching intermittent self catheterisation.



Curriculum References

Undergraduate Medicine

Therapeutic Procedure no. 22, p.14 - Male and female urinary catheterisation - Passing a tube into the urinary bladder to permit drainage of urine, in male and female patients.

Source: GMC Outcomes for Graduates: Tomorrow's Doctors

Postgraduate Nursing

Practise catheterisation on model/manikin under supervision until confident.

Catheter care RCN guidance for nurses (RCN) 2012



For demonstration or to order

Limbs & Things Ltd. Sussex Street, St Phillips, Bristol, BS2 0RA, UK

T: 00 44 (0)117 311 0500 E: gbsales@limbsandthings.com

limbsandthings.com